



EFC-2600

User Manual



Reading through accompanying documents is a mandatory action before using this equipment



810-C1101-533 Ver. A

Content

1	Introduction	2
1.1	Outline	2
1.2	Intended Use	2
1.3	Indication for use	2
1.4	Proper Instrument Use	3
2	Safety Information	4
2.1	Displays for Safety Use	4
2.2	Symbols and Labels	8
2.3	Protective Packing Symbols	10
2.4	Product Labels	19
2.5	Service Life	19
3	Instrument description	20
3.1	Introduction	20
3.2	Standard accessories	22
4	Introduction of User Interface	23
4.1	Patient Management	23
4.2	Capture Image	30
5	Specifications	51
5.1	Product Specification	51
5.2	Environmental conditions	53
5.3	Electric Rating	53
5.4	Fuse specification:	54
6	Maintenance	55
6.1	Lens Cleaning	55
6.2	Chinrest and forehead rest	55
6.3	LCD Monitor/Touch panel	55
7	Installation	56
7.1	Forehead adaptor installation	56
8	Software Update	57

1 Introduction

1.1 Outline

The Ezer non-mydratic fundus camera EFC-2600 is an auto 3D tracking, fast and easy to use retinal imaging system. EFC-2600 is designed to provide images of the eye as an aid to clinicians in the diagnosis of diabetic retinopathy, AMD, glaucoma and other retinal diseases.

1.2 Intended Use

EFC-2600 provides non-mydratic color retina and external images of the eye as an aid to clinicians in the evaluation and diagnosis of eye disease. EFC-2600, which uses NIR LED as illumination during alignment to the retina of patients' eyes, provides images only and does not provide any diagnostic or pathological analysis or conclusions in any way. EFC-2600 is also not intended to be used as the sole diagnostic aid in disease identification, classification or management.

1.3 Indication for use

EFC-2600 is a non-contact, high resolution digital imaging device which is suitable for photographing, displaying and storing images of the retina and external areas of the eye to be evaluated under non-mydratic conditions.

EFC-2600 is indicated for in-vivo viewing of the posterior and external area of the eye and the images are intended for use as an aid to clinicians in the evaluation, diagnosis and documentation of ocular health.

1.4 Proper Instrument Use

1. EFC-2600 is a medical device; it must be operated by properly trained and qualified person(s) only. The operation should be supervised by a physician.
2. Please be sure to read the user manual to understand the safety precautions before operating this device.
3. Always enter patient information first.
4. Prepare patient contact surfaces (forehead and chin rest) according to the cleaning method in this manual.
5. Instantly turn off the power switch of this instrument and disconnect the power cable if uncertain problems arise.
6. Clean ocular lens frequently to ensure good image quality.
7. Adjust the height of motorized adjustable table properly to ensure patient's comfort during the examination.
8. Align the patient's eye position to the canthus indicator mark on the chin and forehead rest assembly.
9. Dim the room lights to allow natural dilation of the patient's pupil and to provide a comfortable visualization of the fixation target without glare.
10. Inspection of the systems' functionality should be conducted once a year or whenever any repairs are made.

2 Safety Information

2.1 Displays for Safety Use

<u>Display</u>	<u>Meaning</u>
 WARNING	"WARNING" indicates the presence of a hazard that could result in severe personal injury.
 CAUTION	"CAUTION" indicates the presence of a hazard that could result in minor injury.
NOTE	"NOTE" provides useful information for operation which is important.

 WARNING	Accessory equipment connected to the digital interfaces must be certified according to the respective IEC standards (e.g., IEC 60950 for laptop or IEC 60601-1 for medical equipment). Furthermore, all configurations shall comply with the system standard IEC 60601-1-1 and IEC 60601-1:2005. Any person who connects or installs devices to the system has responsibility to verify that compliance. If in doubt, consult the Ezer local representative or distributor.
 WARNING	To avoid risk of electric shock, this equipment must only be connected to the supply mains with protective earth.
 WARNING	Do not modify this equipment without authorization of the manufacturer.
 WARNING	THE EFC-2600 CANNOT REPLACE CLINICAL JUDGEMENT AND IS INTENDED TO BE USED ONLY IN CONJUNCTION WITH OTHER CLINICAL TOOLS CONSIDERED TO BE THE

	STANDARD OF CARE FOR MEASUREMENT AND DIAGNOSIS OF THE EYE.
 WARNING	The EFC-2600 is a medical device. The software and hardware has been designed in accordance with U.S., European and other international medical device design and manufacturing standards. Unauthorized modification of the EFC-2600 software or hardware, or any addition or deletion of any application in any way can jeopardize the safety of operators and patients, the performance of the instrument, and the integrity of patient data. <u>Any changes, additions or deletions to factory installed applications, operating system or modifications to hardware in any manner VOIDS the Warranty completely.</u>
 WARNING	EFC-2600 is not intended for home use and may not be stored or operated in environment conditions other than those prescribed. (see <i>Specification</i>)
 WARNING	Phototoxicity Because prolonged intense light exposure can damage the retina, the use of the device for ocular examination should not be unnecessarily prolonged, and the brightness setting should not exceed what is needed to provide clear

	visualization of the target structures. The retinal exposure dose for a photochemical hazard is a product of the radiance and the exposure time. If the value of radiance were reduced in half, twice the time would be needed to reach the maximum exposure limit.
 WARNING	Do not obstruct the mains power switch or position the equipment where the connection to the mains line can be accidentally disconnected.
 WARNING	Equipment is not suitable for use in the presence of a Flammable Anesthetic Mixture with Air, Oxygen, or Nitrous Oxide.
 WARNING	The EFC-2600 has no special protection against harmful ingress of water or other liquids (classified IPX0). To avoid damage to the instrument and cause a safety hazard, the cleaning solutions, including water, should not be directly applied to the device. Using a dampened cloth (without dripping), is a good method to clean the exterior surface of the enclosure.
 WARNING	The patient cannot touch any electrical device that is not powered by EFC-2600 with any part of his or her body while being examined. In addition, the EFC-2600 operator must not attempt to touch the patient and any electrical device that is not powered by EFC-2600 at the same time while examining the patient. Failure to do so could result in electrical shock to the

	patient and/or operator.
 WARNING	Do not connect the instrument with anything other than specified. Otherwise, it may result in fire or electric shock. For details of purchasing accessories, please contact a Ezer representative or distributor.
 CAUTION	Be sure to hold the bottom of the base when EFC-2600 is moved.
 CAUTION	Do not operate the LCD monitor with wet hands or hard objects. The LCD monitor may be damaged.
 CAUTION	When adjusting chinrest, be careful not to pinch the patient's hand.
 CAUTION	The device need to install on the stable table. Do not install in location that are unstable or exposed to vibration.
 CAUTION	Federal law restricts this device to sale by or on the order of a Physician or Practitioner (CFR 801.109(b)(1)).

2.2 Symbols and Labels



Presence of electrical shock hazard.

Note: Indicates risk of electrical shock due to the presence of uninsulated high voltage inside the instrument. Do not remove the instrument cover or parts.



Circuit Breaker



Type B applied parts.

Note: This instrument complies with the specified requirements to provide protection against electrical shock, particularly regarding allowable patient leakage current.



Authorized Representative in the European Union

Medical Device Safety Service (MDSS) GmbH

Schiffgraben 41, 30175 Hannover, Germany



Non-mydriatic Fundus Camera (EFC-2600)
User Manual



Serial number



Catalog number / part number



Prescription Use

2.3 Protective Packing Symbols

The protective packing symbols specify the handling requirements and the transport and storage conditions.



Fragile, Handle with care




Keep Dry



This end up



Relative Humidity

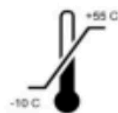

Temperature



DO NOT
STACK

Do not stack



2 Layers only



Waste Electrical and Electronic Equipment (WEEE) Recycling

Instructions

When determined that the device is ready for disposal, it is to be recycled following the policies and procedures reflecting respective country's requirements. Do not dispose of device as general waste.

Product Compliance



Indicates this equipment contains Type B applied parts

The EFC-2600 is classified as follows:



Class I Equipment – Protection against electrical shock.

Type B – Degree of protection against electric shock of applied part (chin and forehead rests).

Ordinary Equipment (IPX0) – Degree of protection against ingress of liquids (none).

Continuous Operation – Mode of operation

Electromagnetic Compatibility (EMC): EN 60601-1-2

The EFC-2600 device has been tested to comply with the emission and Immunity requirements of EN60601-1-2.

The EFC-2600 is intended for use in an electromagnetic environment where radiated RF disturbances are not beyond the standard defined in EN60601-1-2.

Certification: under IEC 60601-1

Guidance and manufacturer's declaration – electromagnetic emissions		
The EFC-2600 is intended for use in the electromagnetic environment specified below. The customer or the user of the EFC-2600 should assure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	The EFC-2600 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class A	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	The EFC-2600 is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.

Guidance and manufacturer's declaration – electromagnetic immunity			
The EFC-2600 is intended for use in the electromagnetic environment specified below. The customer or the user of the EFC-2600 should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) IEC 61000-4-2	6 kV contact 8 kV air	6 kV contact 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/burst IEC 61000-4-4	2 kV for power supply lines 1 kV for input/output lines	2 kV for power supply lines 1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	1 kV line(s) to line(s) 2 kV line(s) to earth	1 kV line(s) to line(s) 2 kV line(s) to earth	Mains power quality should be that of a typical commercial or hospital environment.
interruptions and voltage variations on power supply input lines	<5 % <i>UT</i> (>95 % dip in <i>UT</i>) for 0,5 cycle	<5 % <i>UT</i> (>95 % dip in <i>UT</i>) for 0,5 cycle	Mains power quality should be that of a typical commercial or hospital environment. If the user of the EFC-2600 requires

IEC 61000-4-11	40 % <i>UT</i> (60 % dip in <i>UT</i>) for 5 cycles 70 % <i>UT</i> (30 % dip in <i>UT</i>) for 25 cycles <5 % <i>UT</i> (>95 % dip in <i>UT</i>) for 5 sec	40 % <i>UT</i> (60 % dip in <i>UT</i>) for 5 cycles 70 % <i>UT</i> (30 % dip in <i>UT</i>) for 25 cycles <5 % <i>UT</i> (>95 % dip in <i>UT</i>) for 5 sec	continued operation during power mains interruptions, it is recommended that the EFC-2600 be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE <i>UT</i> is the a.c. mains voltage prior to application of the test level.			


Guidance and manufacturer's declaration – electromagnetic immunity			
The EFC-2600 is intended for use in the electromagnetic environment specified below. The customer or the user of the EFC-2600 should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 Vrms	Portable and mobile RF communications equipment should be used no closer to any part of the EFC-2600, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 1,2 \sqrt{P}$
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2,5 GHz	3 V/m	$d = 1,2 \sqrt{P}$ 80 MHz to 800 MHz $d = 2,3 \sqrt{P}$ 800 MHz to 2,5 GHz where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).

			Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, a should be less than the compliance level in each frequency range.b Interference may occur in the vicinity of equipment marked with the following symbol: 
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures objects and people.			
Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the EFC-2600 is used exceeds the applicable RF compliance level above, the EFC-2600 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the EFC-2600. Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.			

Recommended separation distances between portable and mobile RF communications equipment and the EFC-2600			
The EFC-2600 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the EFC-2600 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the EFC-2600 as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1,2 \sqrt{P}$	80 MHz to 800 MHz $d = 1,2 \sqrt{P}$	800 MHz to 2,5 GHz $d = 2,3 \sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer. NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

2.4 Product Labels

EFC-2600 system labels (Samples Only):

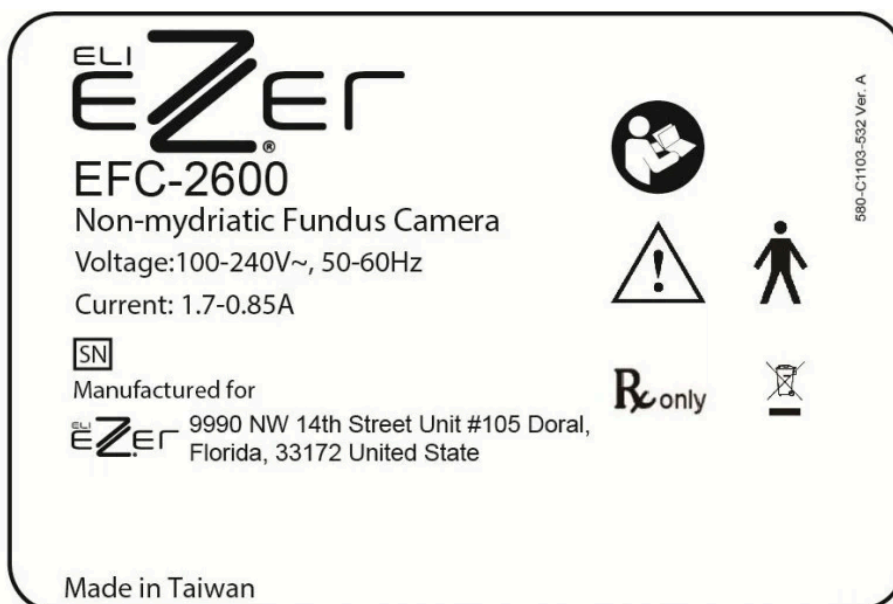


WARNING: Do not connect the instrument with anything other than specified. Otherwise, it may result in fire or electric shock. For details of purchasing accessories, please contact a Ezer representative or distributor.


2.5 Service Life

The service life of EFC-2600 is five years if specified inspections and maintenance are done.

3 Instrument description

3.1 Introduction

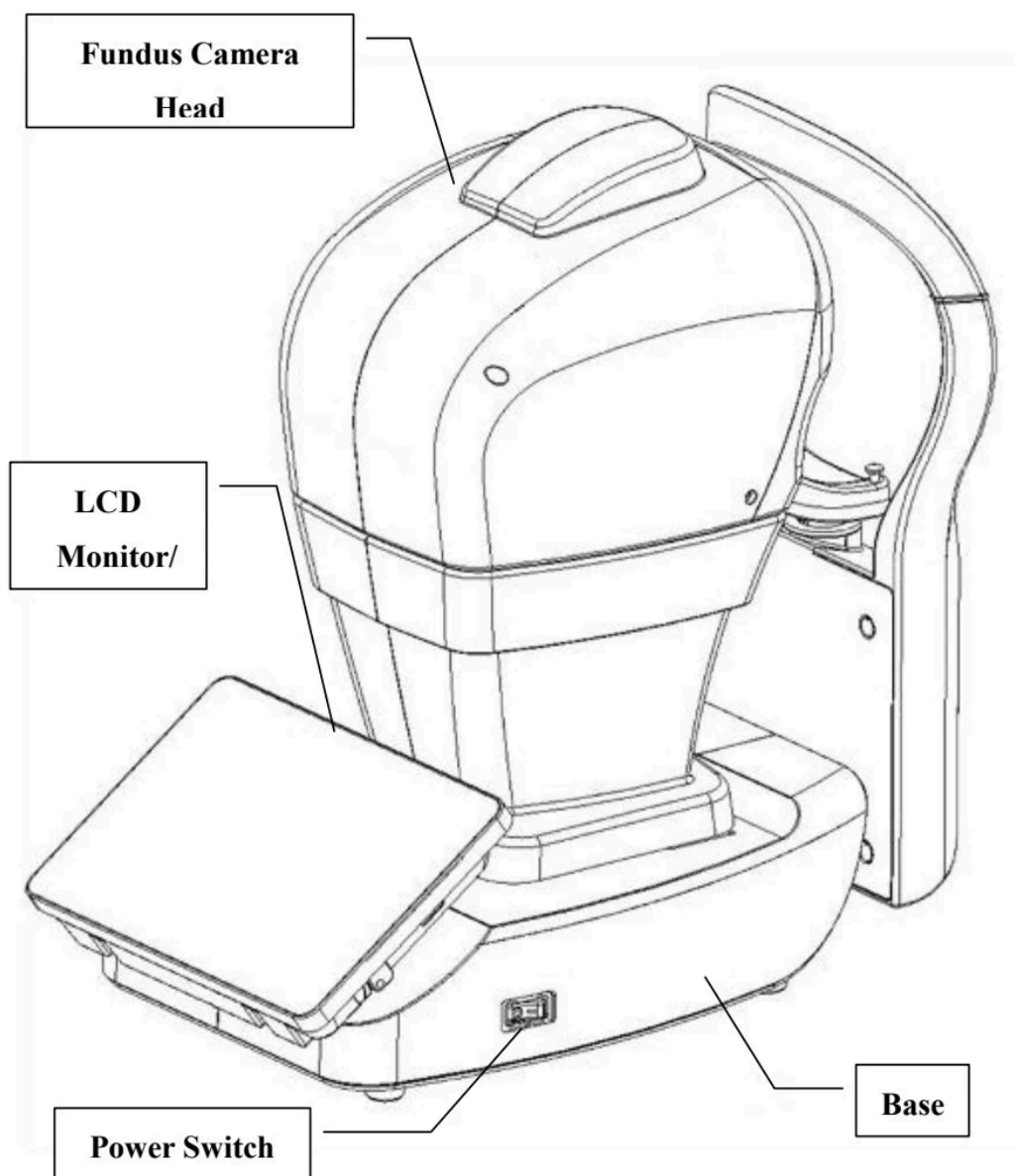

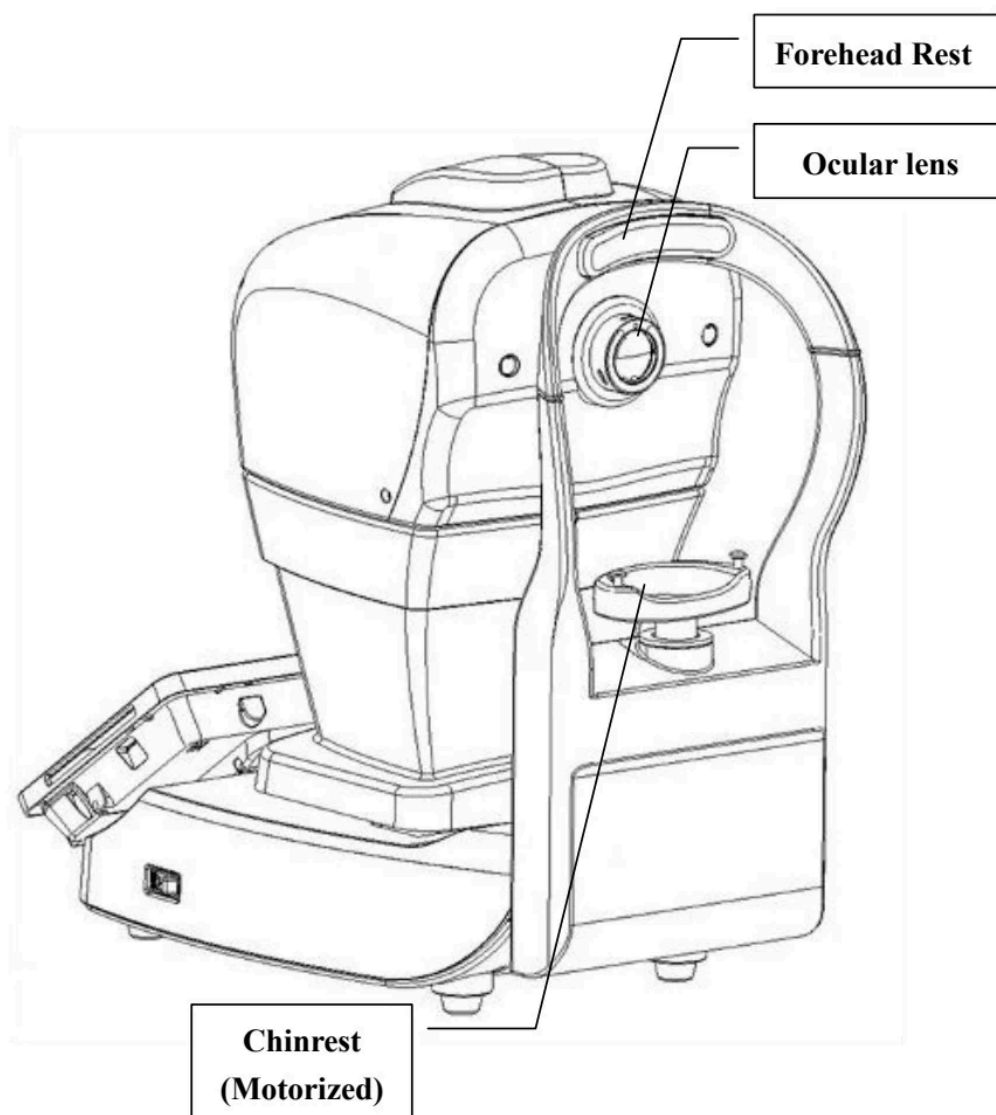




3.2 Standard accessories

Item	Description	Q'ty
1	User Manual	1 pc.
2	AC Power Cord	1 pc.
3	Dust Cover	1 pc.
4	Forehead adaptor	1 pc


4 Introduction of User Interface

Three tabs are displayed on screen after power on and describe in below sections.

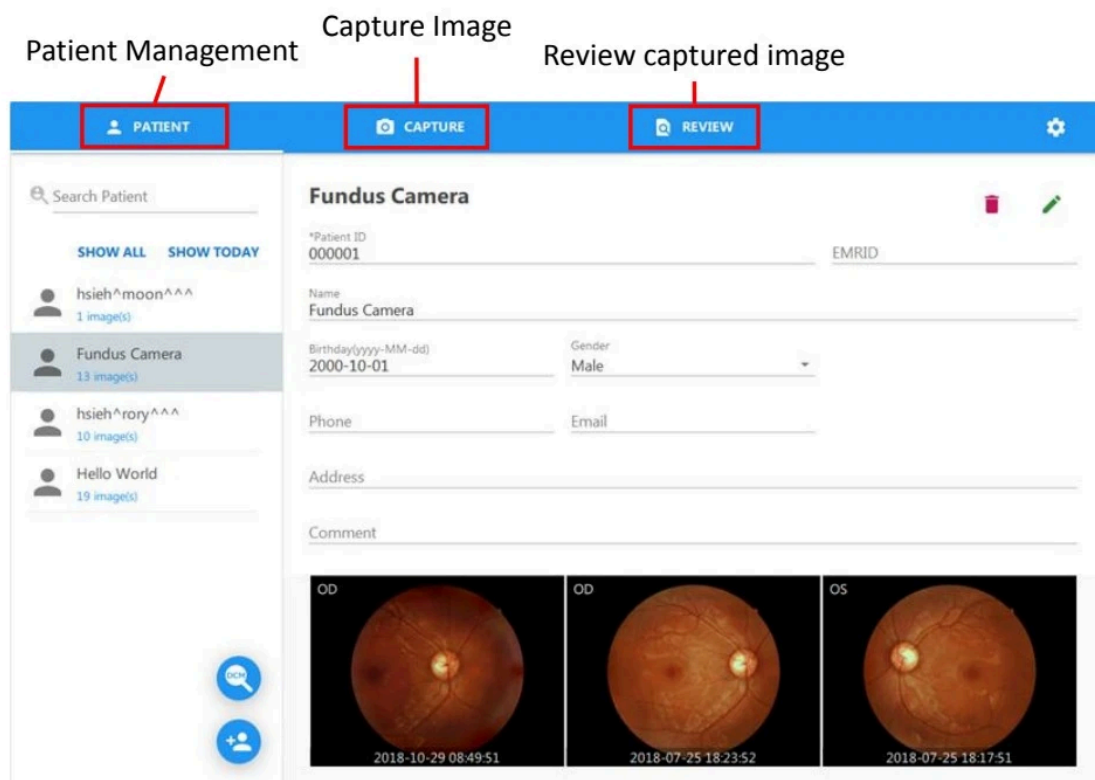


4.1 Patient Management

Select patient, add patient and modify patient information

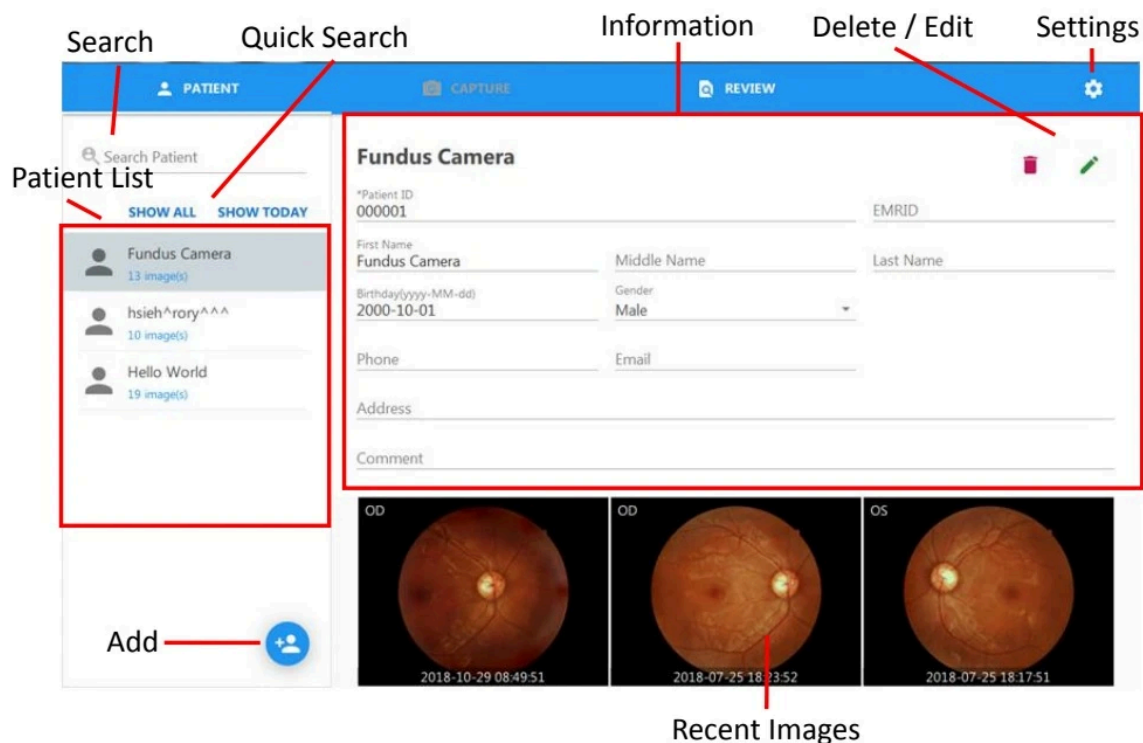


- Patient List: Displays all the patients associated with the search result.

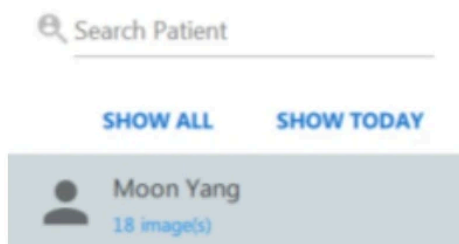


- Search: Provides patient search function by entering keywords.
 - Search by all fields: Fill in keyword such as "Moon"
- Add Patient  : Click to add new patient profile and information.




- The columns with * mark are mandatory before adding new patient to database.
- By clicking , the new patient information is saved to the database. The main screen will then return to the main Patient Information window and the newly added patient is listed under the patient list window. Select the new patient added and patient detail will be displayed accordingly.
- By clicking , it will return to main patient window without saving.
- Patient Information
 - Displaying selected patient information
 - Edit Patient Information  : Click for patient information editing and comment editing saving. The operations are just like Add Patient.
 - Delete Patient  : Click for delete patient information and images. When the “Delete” dialog appears, click “OK” to delete patient information or click “Cancel” to quit.

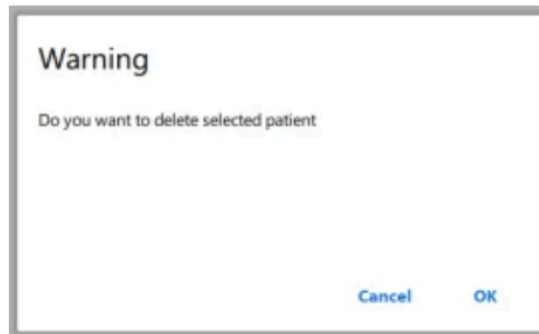


- Settings  : Click to show settings of EFC-2600 program
- General: Enabling Use Cropped Image, the resolution of captured image will be saved 2800x2800

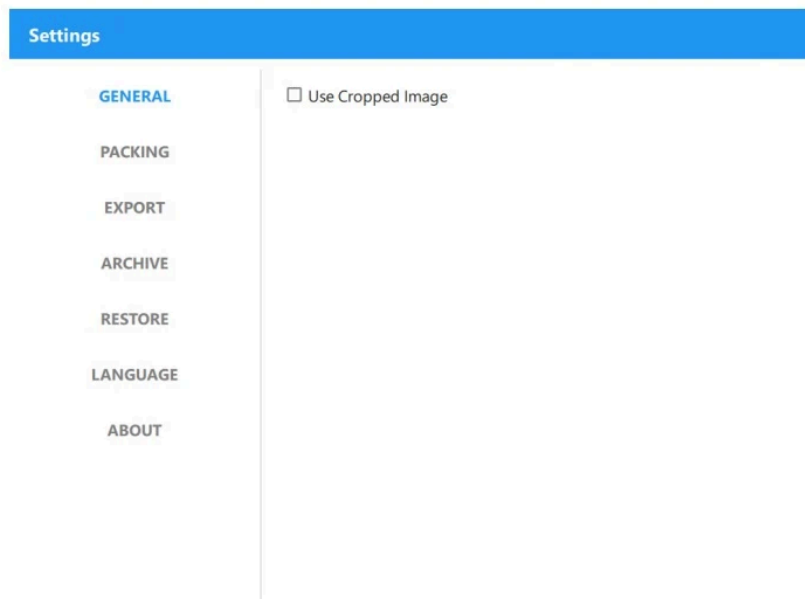


- Packing: Enabling Auto Packing, the camera head moves to packing position before EFC-2600 shutting down.

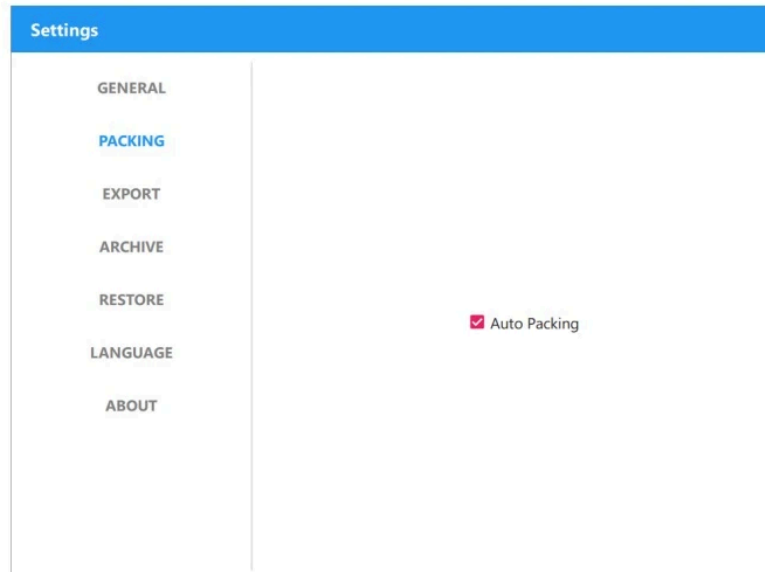


- Export: Enabling Auto Export, the captured image will be copied to specific path automatically where the user set

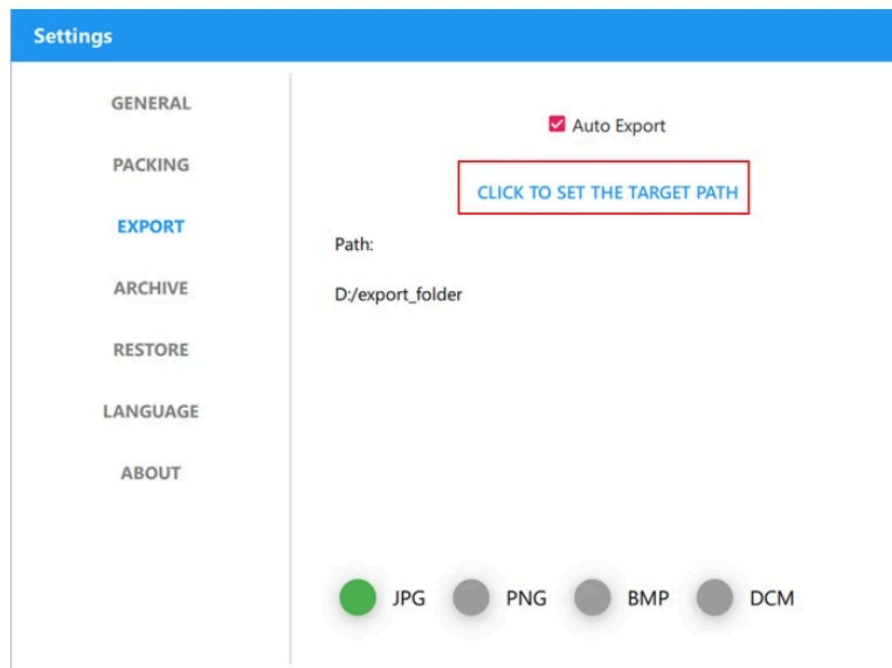

- Archive: Creates a backup file of EFC-2600 database.
Step1. Choose the target folder for backup file

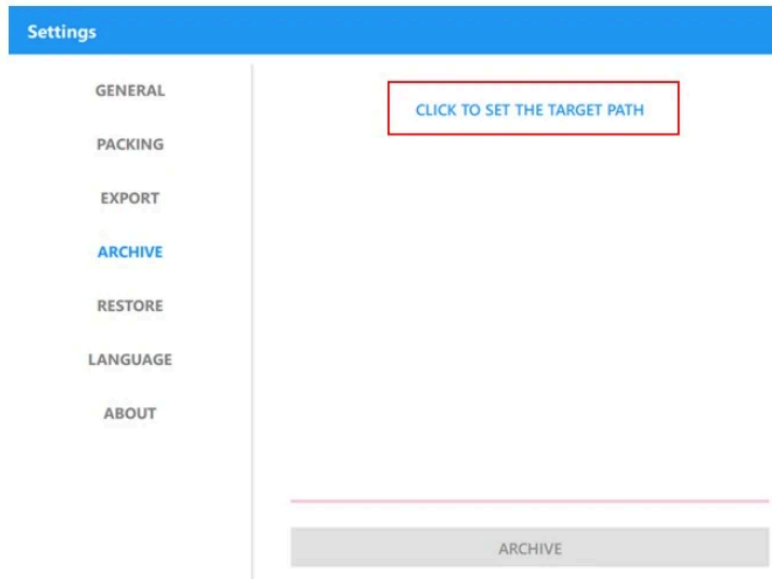


Step2. Click “ARCHIVE” button to create backup file

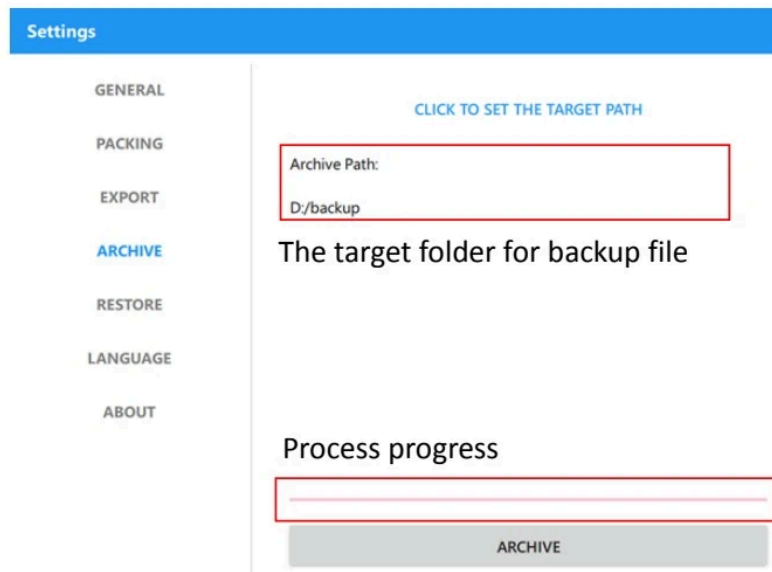


NOTE: The file system of target drive should be NTFS

- Restore: Restore the database of EFC-2600 by backup file


Step1. Choose the target folder where the backup file is saved

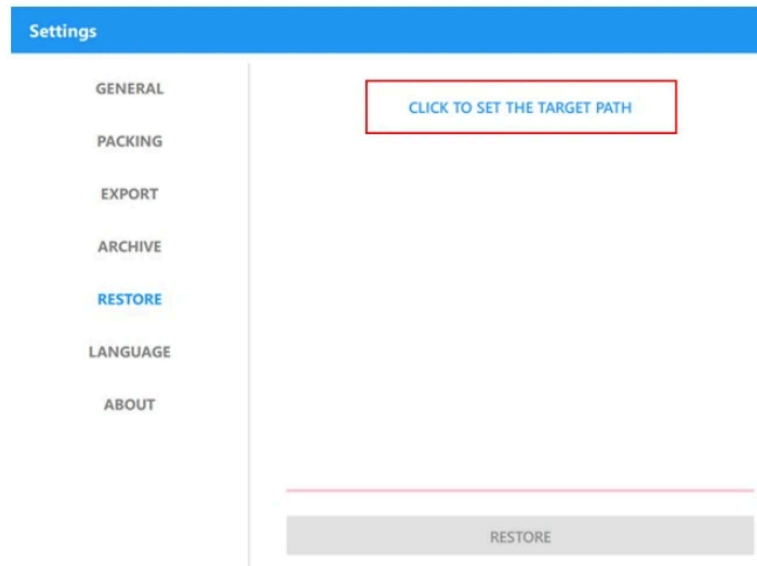


Step2. Choose the backup file from the list and click the “RESTORE” button to restore the database of EFC-2600

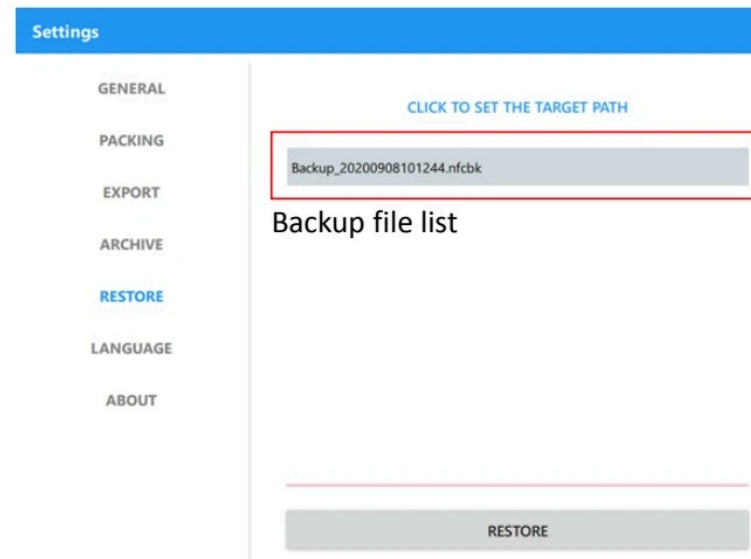

- Language: Multi-labguage selection

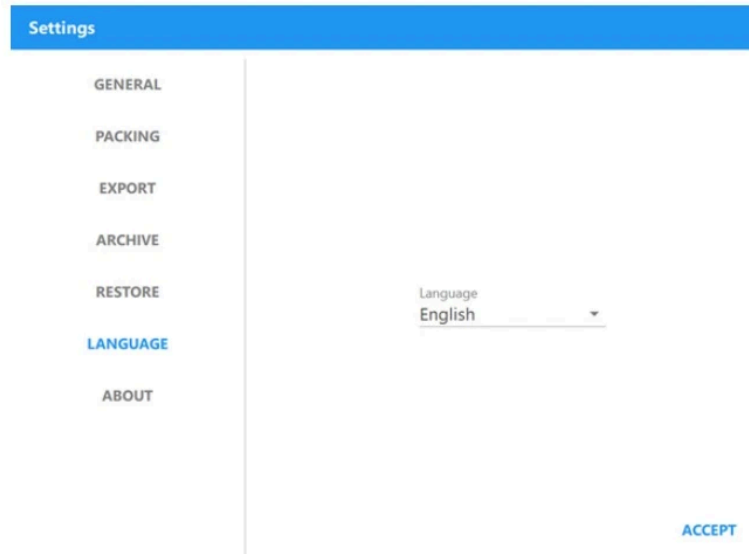


- About: Display software version of EFC-2600

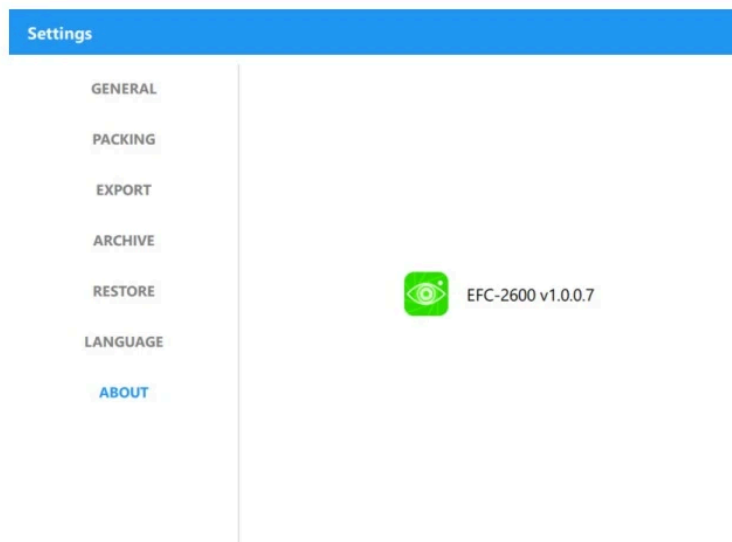

4.2 Capture Image

Main window for image capture

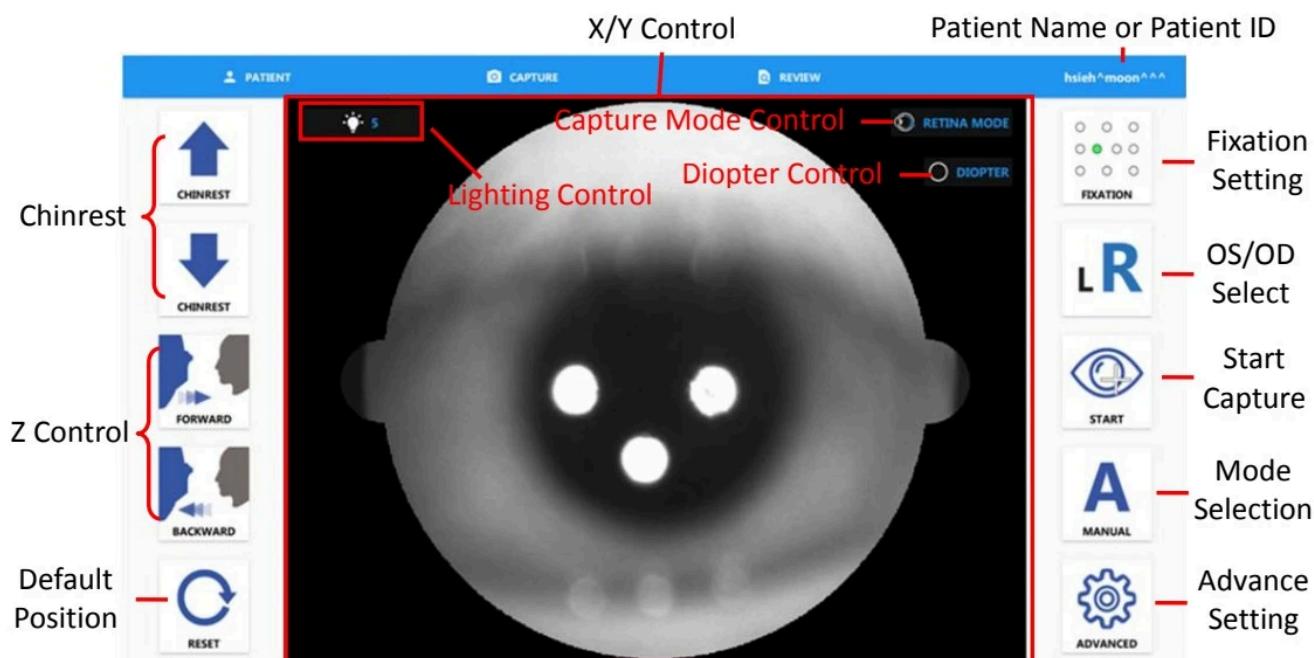

Control Buttons on the screen


- Chinrest: Control Chinrest up and down
- Z Control: Move camera forward or backward
- X/Y Control: Click center of pupil on screen to alignment
- Reset: Reset camera to default position
- Fixation Setting: Selection of the fixation position
- OS/OD Select: Choose OD or OS for image capture
- Start: Click for automatic eye alignment and image capture
- Mode Selection: Auto alignment mode or manual alignment mode
- Advanced: Display settings of EFC-2600 program
 - Semi Auto Mode: User needs to switch eye manually.
 - Full Auto Mode: The EFC-2600 will switch to another eye automatically.
 - Enable manual mode helper: It can help to find suitable working

distance and do capture processes automatically if enable it.


Settings

Mode		☒ Semi Auto	☐ Full Auto
Manual Mode		☒ Enable manual mode helper	


- Capture Mode Control: Display capture mode setting panel and provide retina (default) and cornea mode.

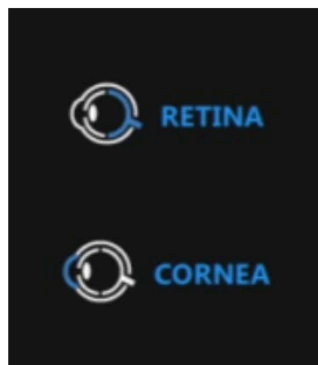

- Diopter Control: Display diopter setting panel and provide empty lens (default), - lens and + lens

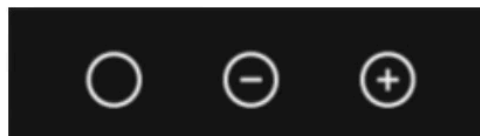

- Lighting Control: Display lighting setting panel and provides viewing and flash LED level

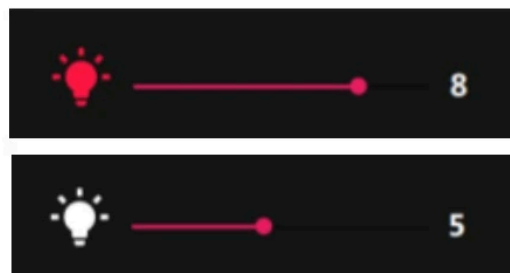


- Dilate Selection: The captured image will save with this setting.

4.2.1 Automatic alignment and focus operation


- Adjust Chinrest with   and table to suitable position.
- Click pupil position on screen to align the camera.
- Click  to start tracking and capture.

4.2.2 Montage mode


- Press and hold the Start button for 3 seconds, the Montage button will be shown.

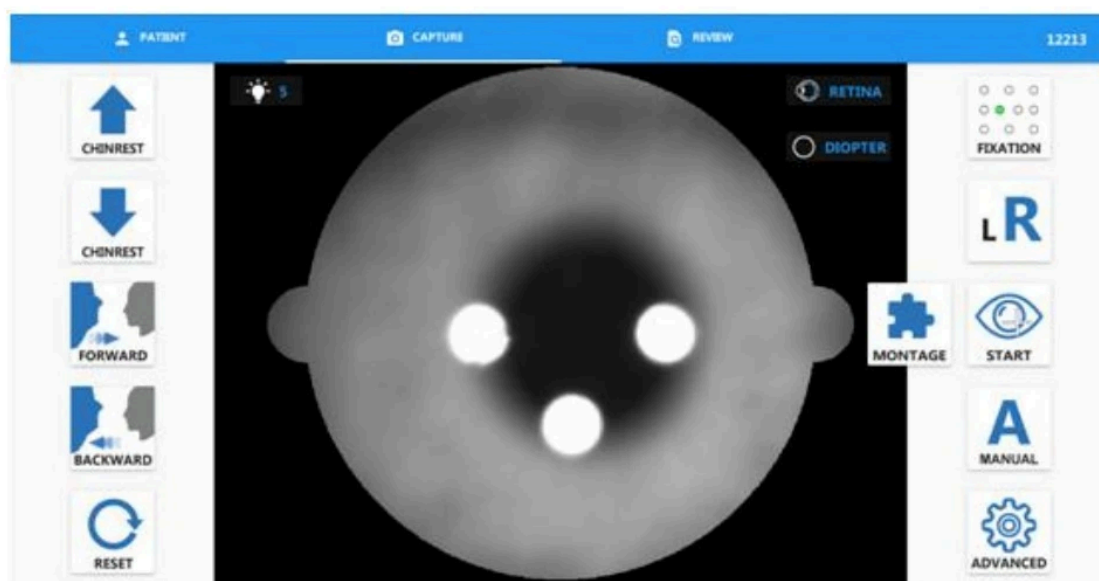


- Click the Montage button and select fixation LEDs.

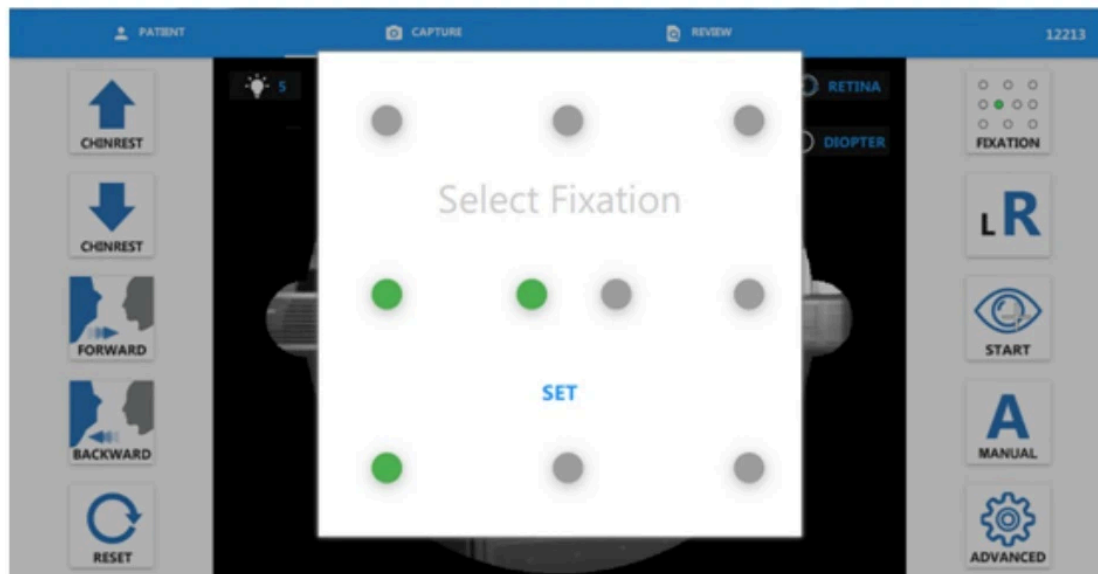

- Click the SET button.
- Click the OK to start montage capture

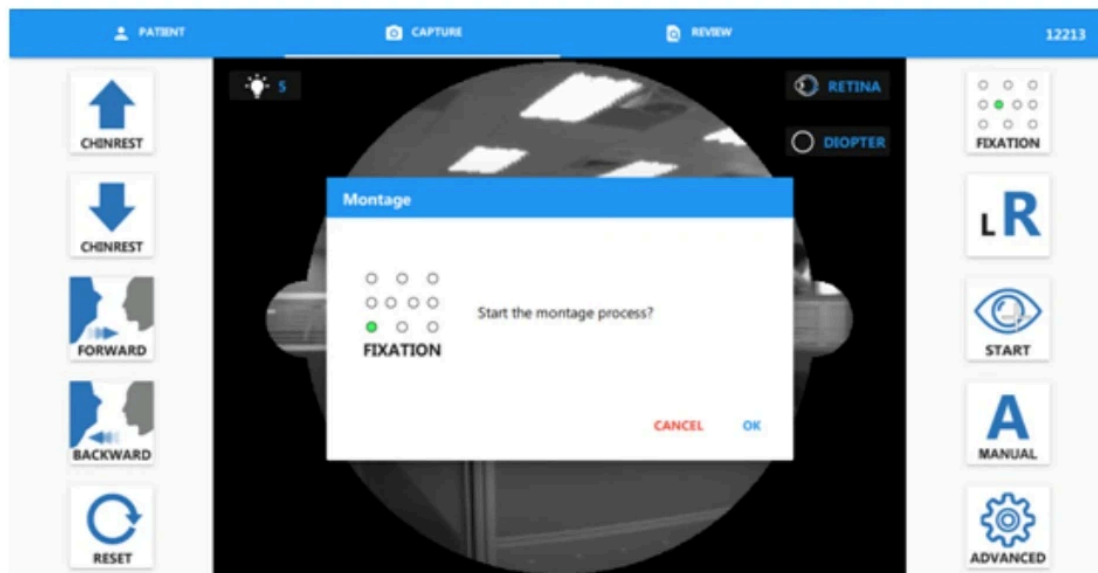

- Click CONTINUE button for next shot.

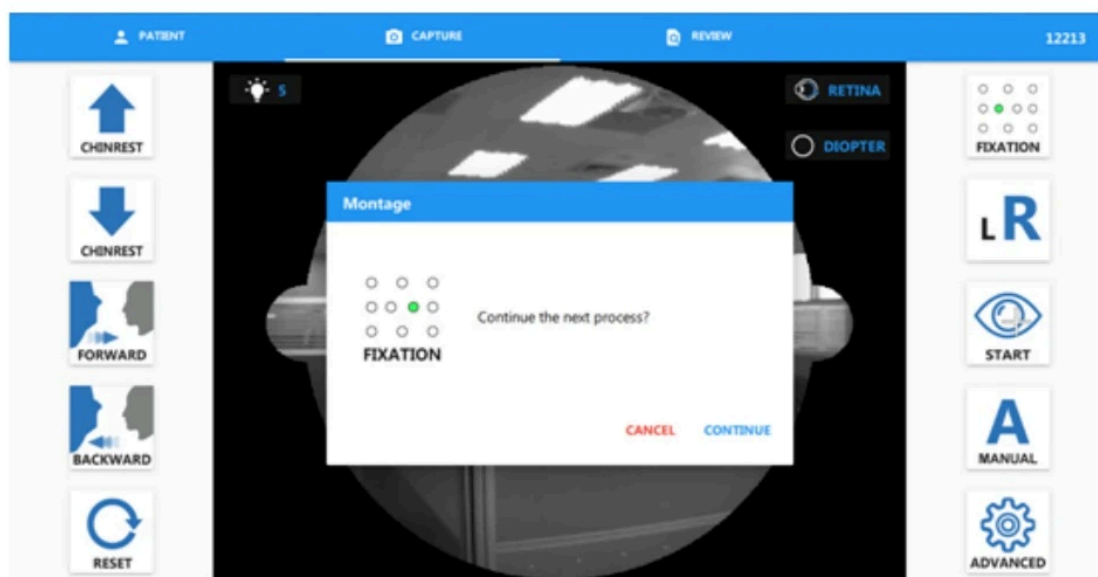


4.2.3 Manual alignment and manual capture


- Click Mode Selection button to manual mode

- Adjust Chinrest by   and table to suitable position.


- Click pupil position on screen or press     buttons to align the camera.

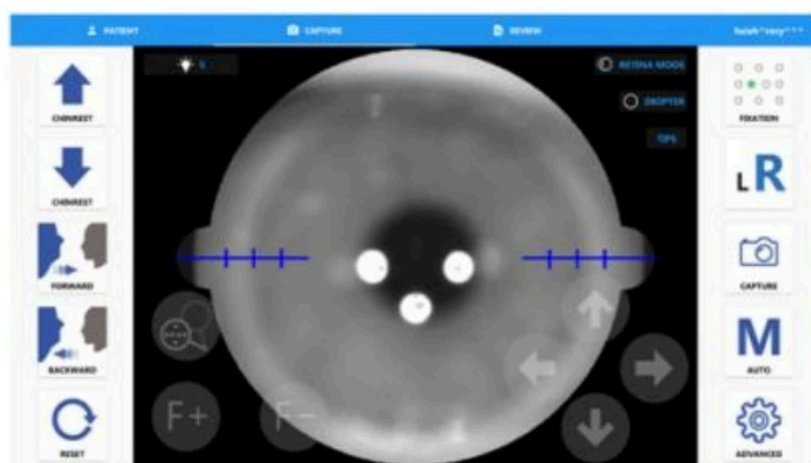


- Approach to the pupil by   buttons until split bar is visible

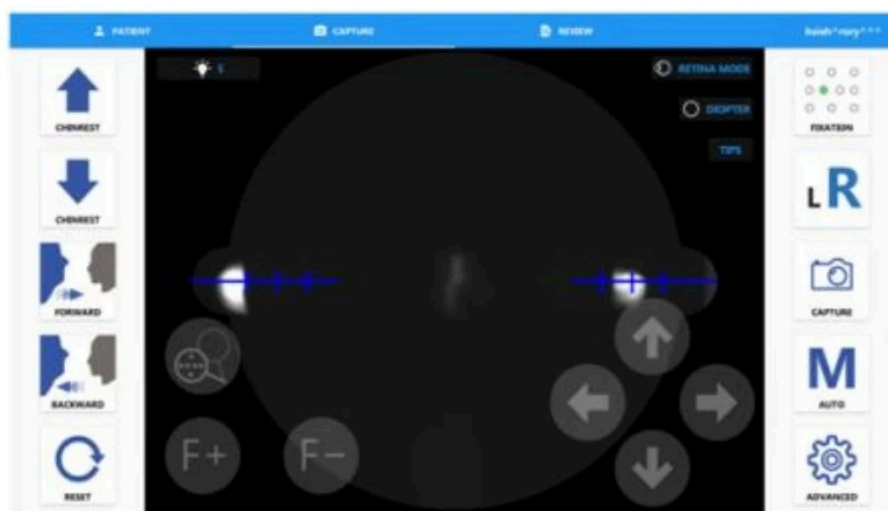


- Align split bar by   (Click on the  button, the fixation mask plate will be removed. User can see the live retina video.)

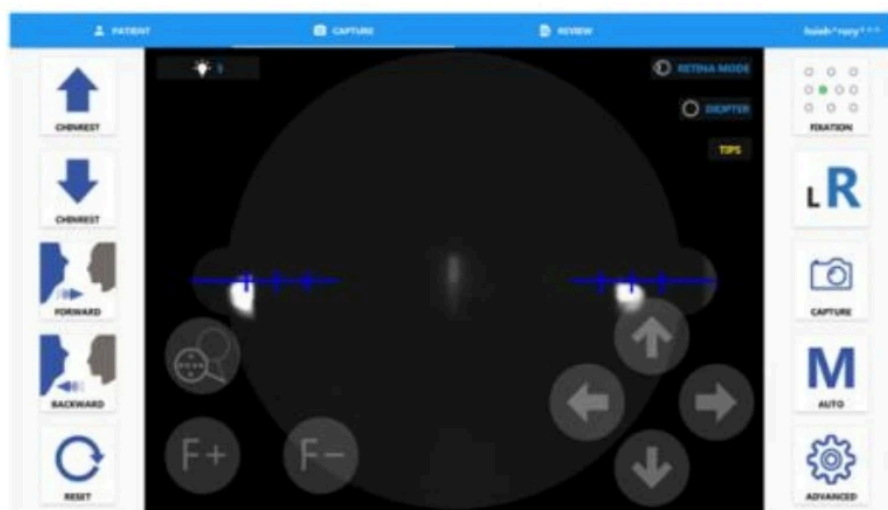

- Enter the pupil by   buttons until two spots are appered

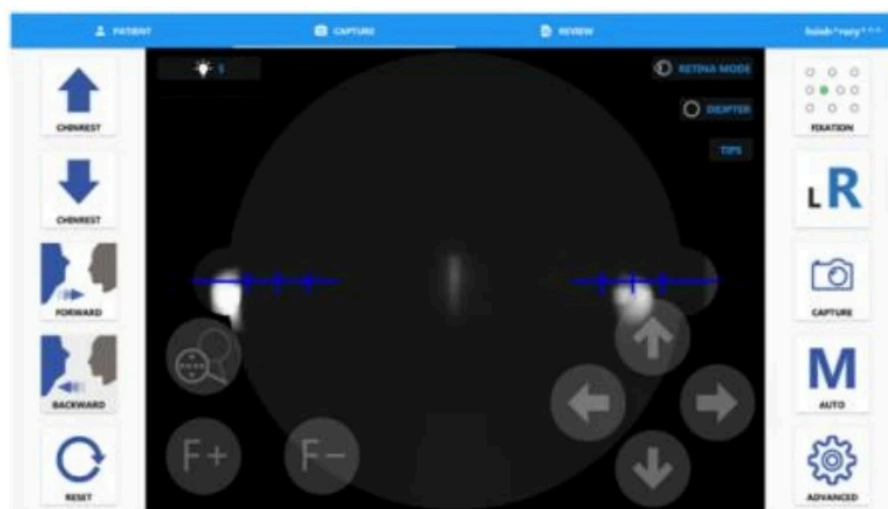


- Use       buttons to align two spots and click button  to capture

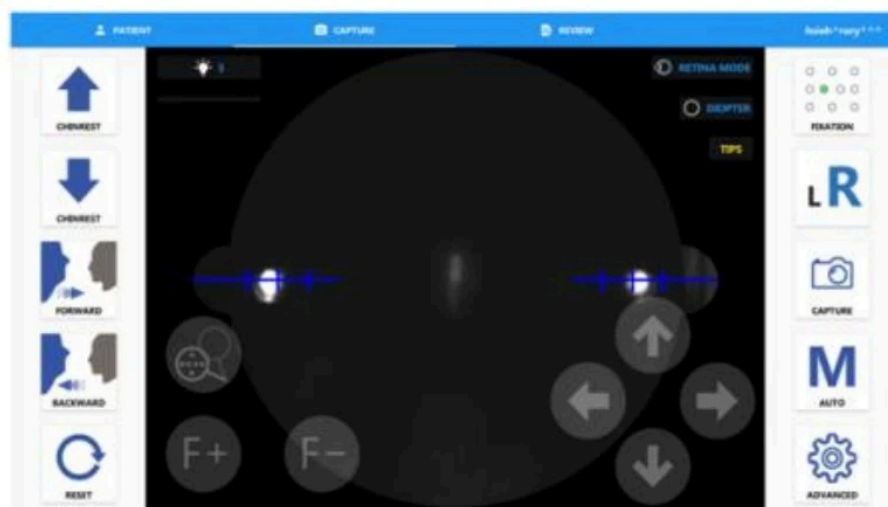

- Click on the **TIPS** button, the operation tips will be shown.


Manual Mode Tips

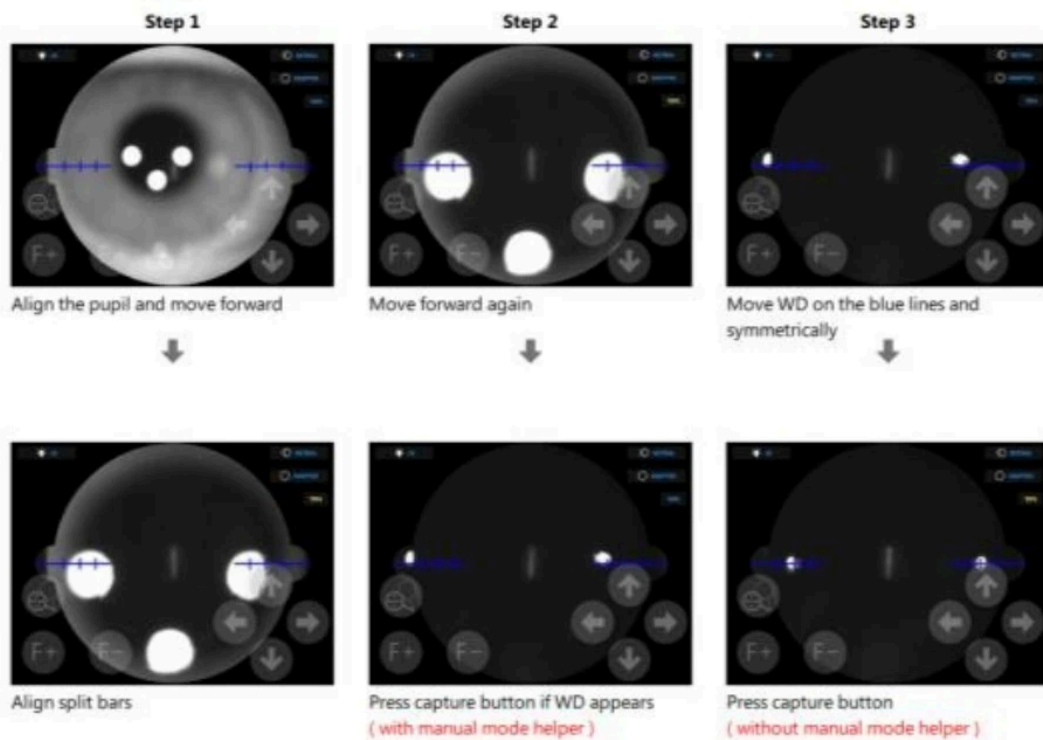

4.2.4 Cornea Capture Mode

- Add forehead adaptor to forehead rest (refer to item 7.1 forehead rest installation)

- Click the  button and then select the cornea capture mode.

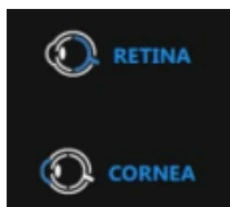

- Adjust Chinrest with   and table to suitable position.

- Click pupil position on screen or press     buttons

to align the camera.




- Use  buttons until the image is clear.



- Click  to capture


4.2.5 Send DICOM Image

When the selected patient which is found from DICOM server, the DICOM sending window is displayed before user leaves the Capture Image page. User can choose images and send them to DICOM server.

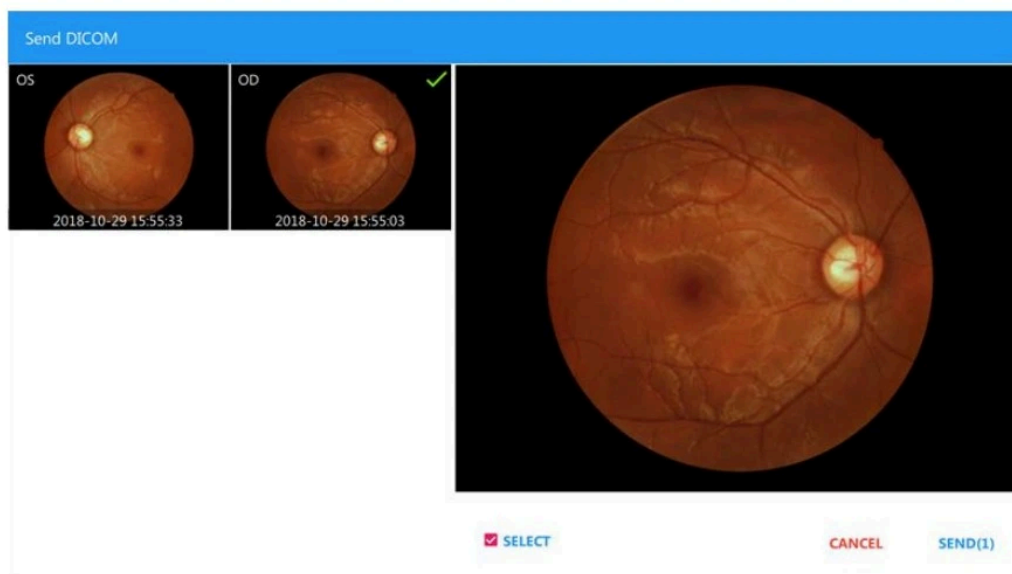


4.2.6 Review Captured Image

Display the Visit List classified by capture date, and the capture time is also displayed with each image. Operator can review images in this window on a particular date.

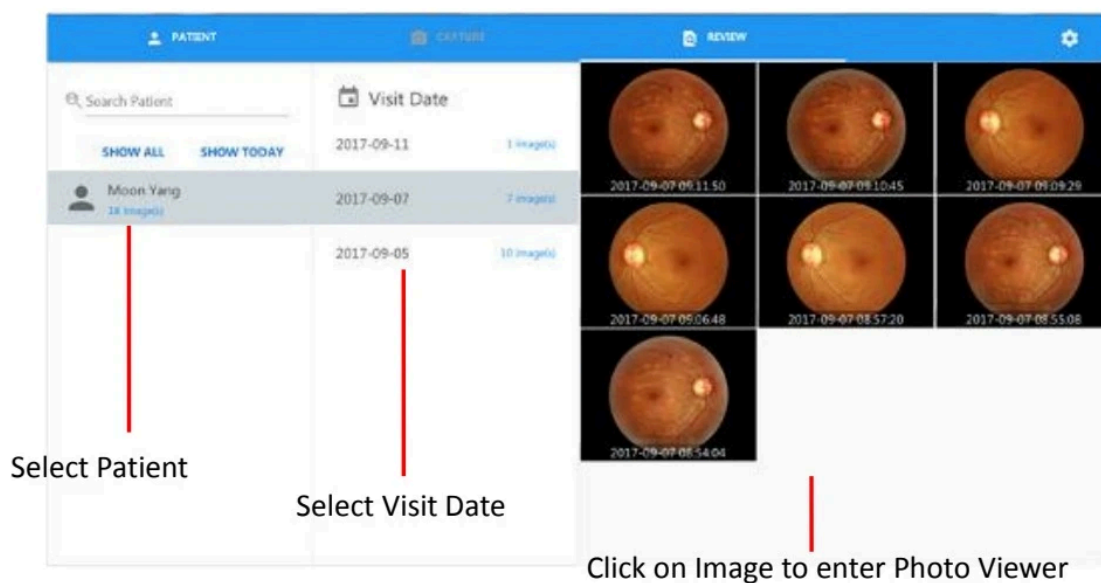


■ Multiple Selections

- Press and hold on image thumbnail can enter the multiple selection mode.
- Click image thumbnail to select image for batch delete or export.

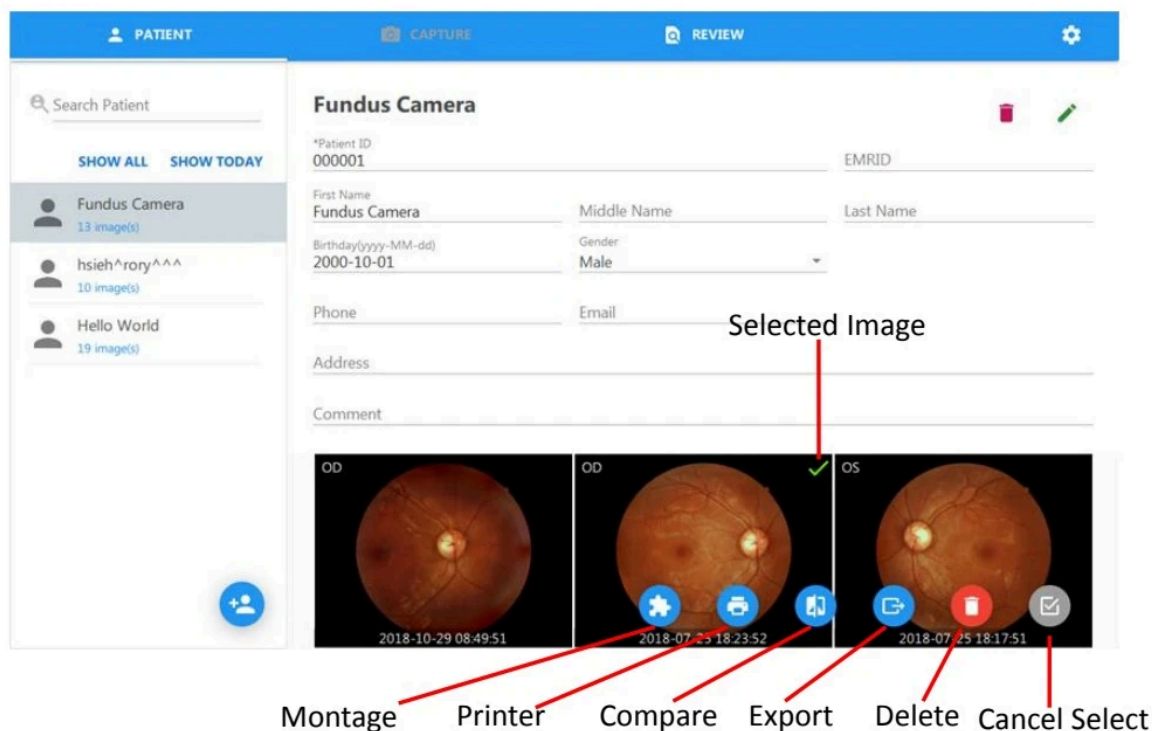

■ Image Delete

- Click on  button, warning dialog will show. Click yes to delete selected image. Becarful, Deleted image cannot restore!

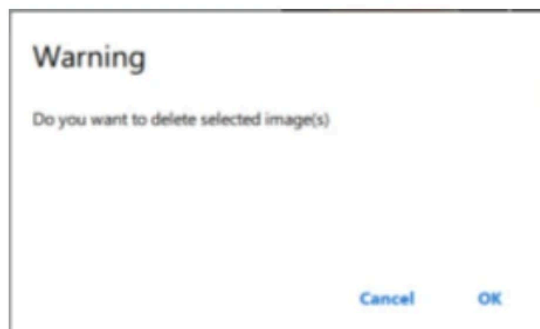


■ Montage


- Click on  button, montage dialog will show.



- Click  button, the visit date list of selected patient will show.

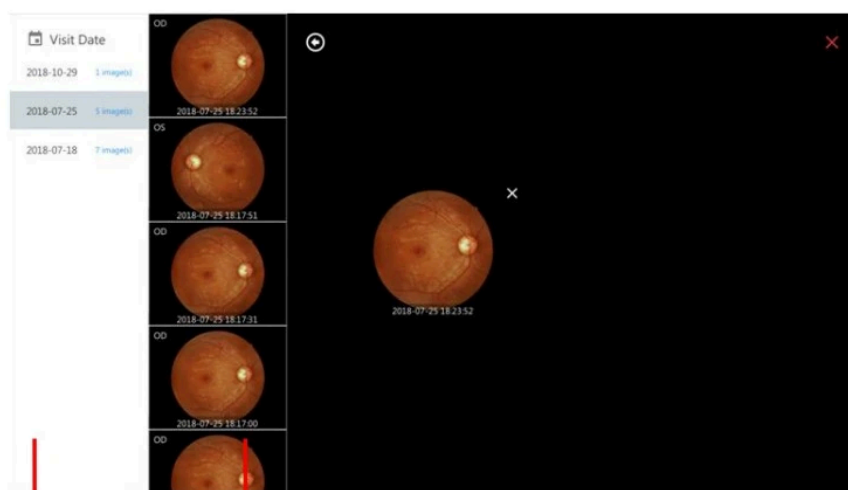


Visit Date List Click on image to assign it to candidate list

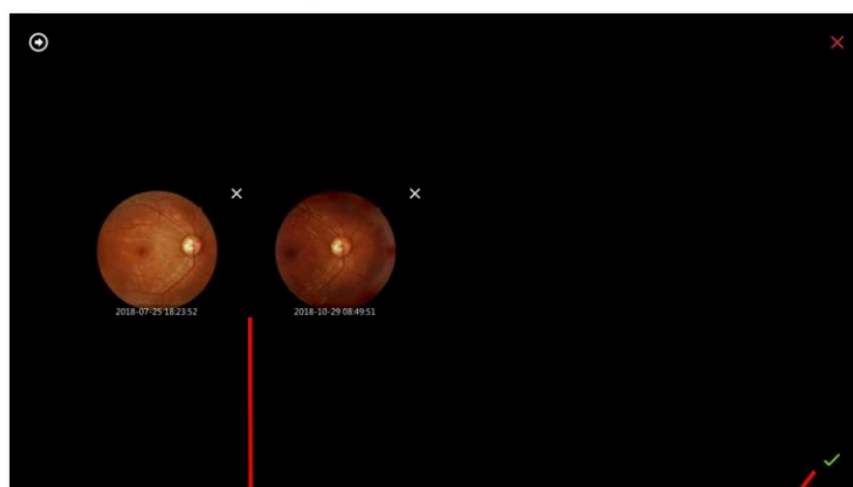

Candidate List


Start Montage

- Click  button and the montage result will be shown.

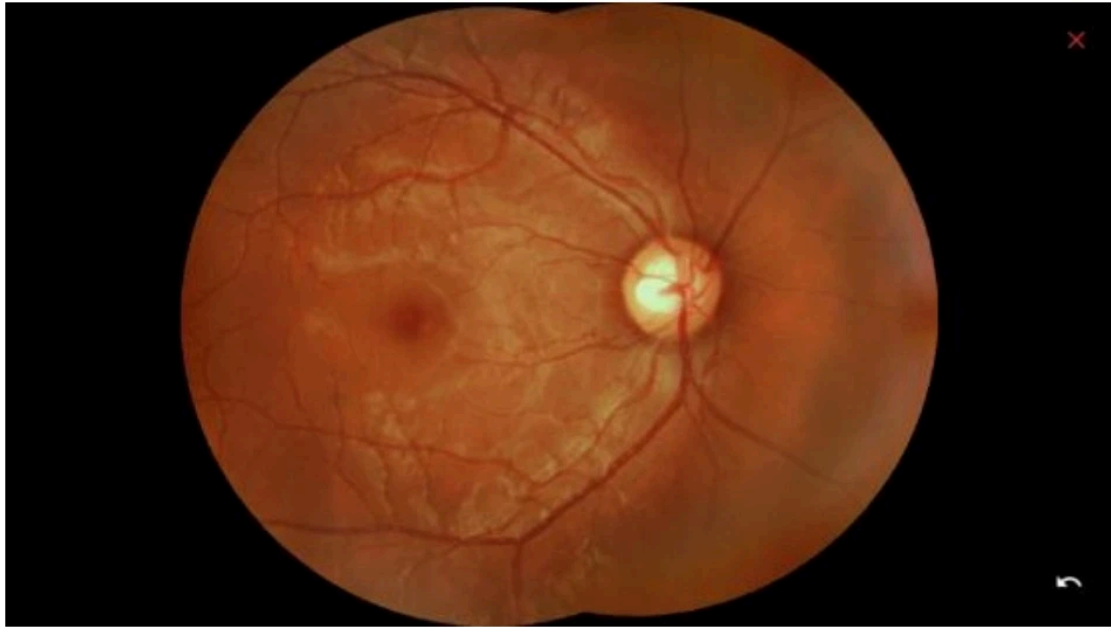

■ Image Export

- Click on  button, warning dialog will be shown and then choose the destination folder, image compression format. The Export ID means the file name should be included patient ID.

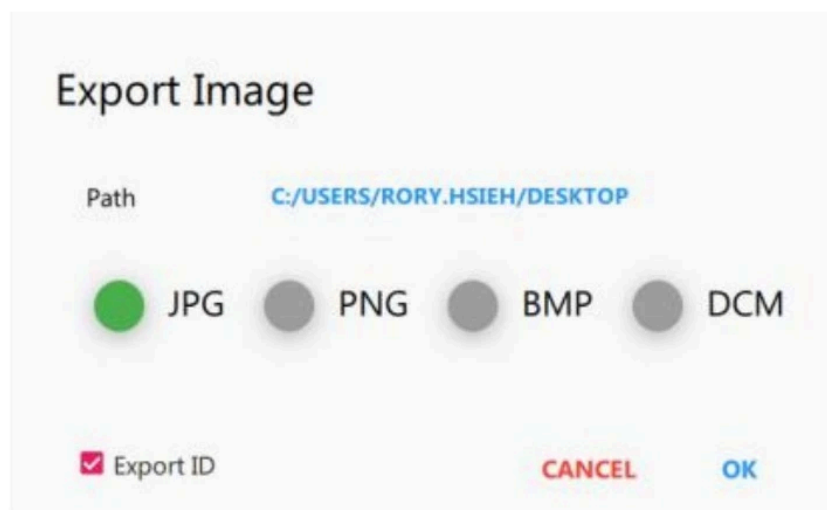


- Press and hold the patient name from patient list, the export button will be appeared. Click the  button and finish export options,


the all images of selected patient will be exported.

■ Image Comparison


- Click on  button, comparison dialog will be shown

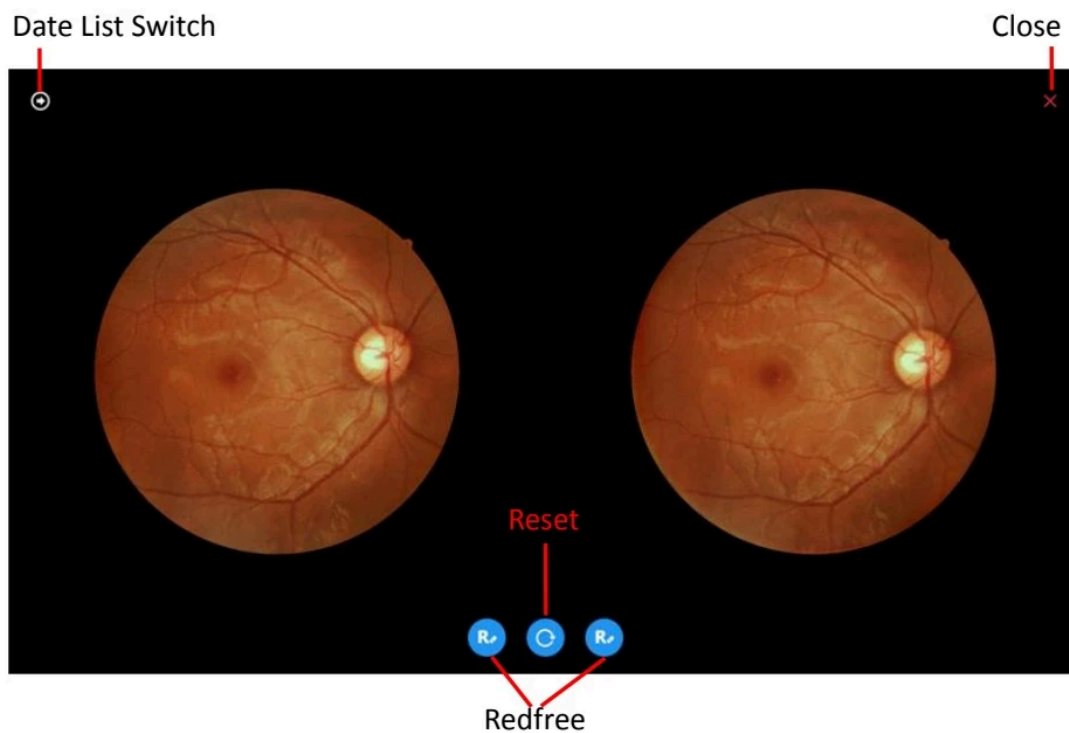

- Date List Switch: Click  button, the visit date list of selected patient will be shown.

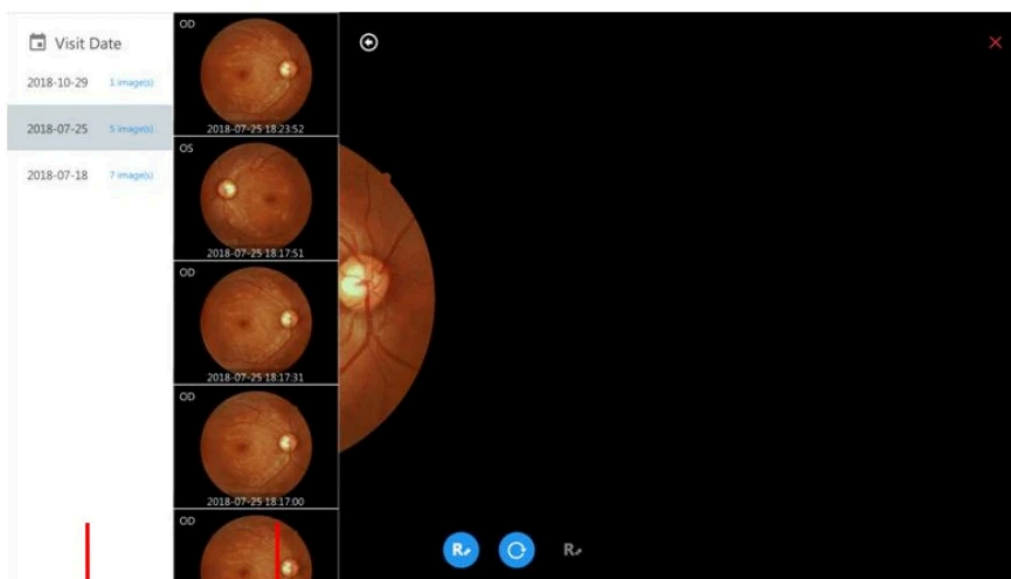


Visit Date List Click on image to assign it to left or right frame

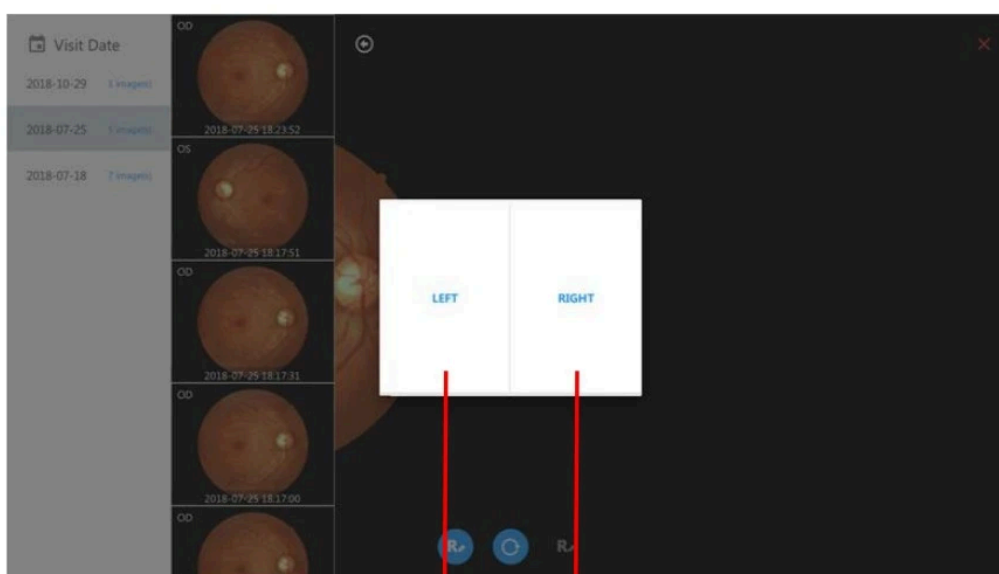


Left frame Right frame

- Reset: Click  button, all image settings are reset to default
- Redfree: Click  button to show the redfree image
- Printer
 - Click on  button, printer dialog will show. Click the PRINT to print selected image.

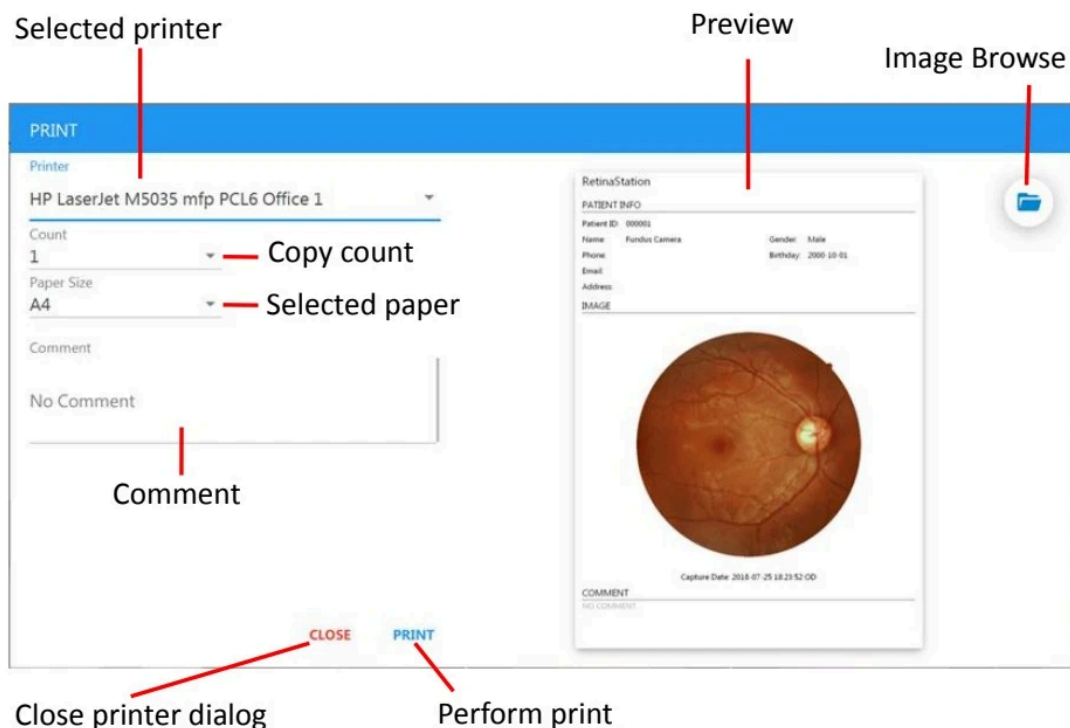


- Select two images and click on  button, printer dialog will show. Click the PRINT to print selected images.

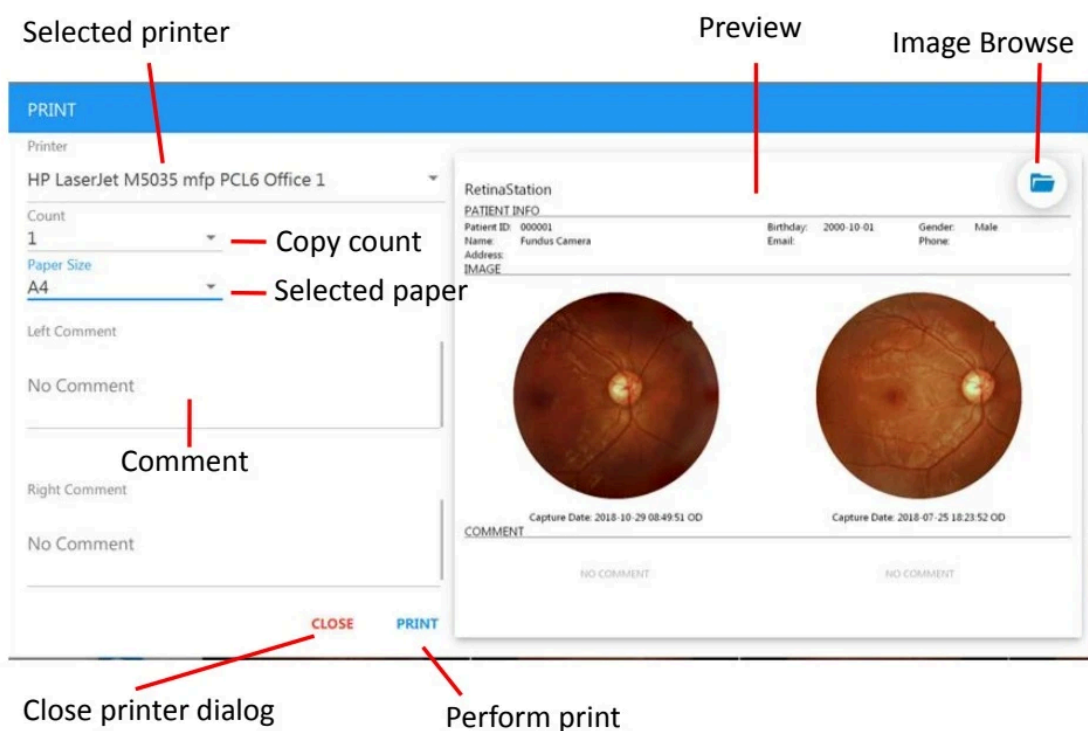

4.2.7 Photo Viewer

- Click image of Image List, the Photo Viewer will show. The Photo Viewer includes below functions.

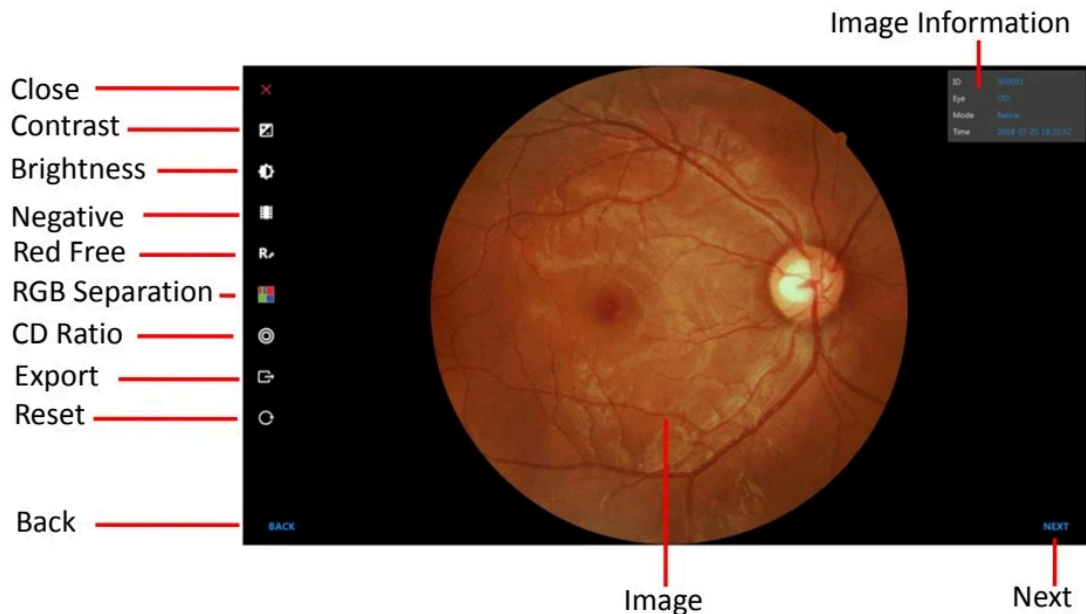


- Image Information
ID: patient ID
Eye: captured eye
Mode: capture mode
Time: capture date and time
- Close: Exit the Photo Viewer
- Contrast: Adjust the selected image. The  button is contrast level up and the  button is contrast level down.



- Brightness: Adjust the selected image. The  button is brightness level up and the  button is brightness level down.




- Negative: To do negative process for selected image
- RedFree: Remove the red channel and convert to grayscale for selected image.
- RGB Separation: To do RGB channels separation and display in color or grayscale mode. The  button is grayscale mode and the  is color mode.

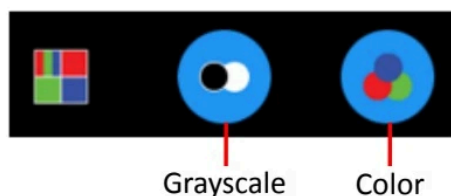


- CD Ratio: To do disc and cup measurement. Below describes how to measure the cup to disc ratio.
Step1. Use the pinch gesture to zoom in/out for region of interest
Step2. Click the CD Ratio button and click the  button.

Step3. Tap the border of disc until the blue line fit it.

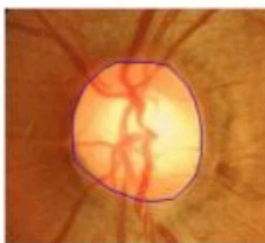

Step4. Click the CD Ratio button and click the  button.

Step5. Tap the border of cup until the yellow line fit it.

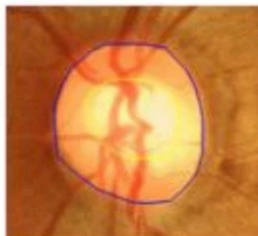

Step6. Click the CD Ratio button and click the  button to save the measurement result. The measurement is displayed on the top-right side of Photo Viewer.

Cup to disc horizontal ratio —
Cup to disc vertical ratio —

ID	000001
Eye	OD
Mode	Retina
Time	2018-07-25 18:23:52
CD H	0.7
CD V	0.6

Click the  button is for clear measurement.

- Export: Export the current image
- Reset: Reset all measurement, scaling ratio, image position and



processing level etc...

Back/Next: Choose the next or previous image

5 Specifications

5.1 Product Specification

Function	Value/Type	Remark
Fundus image		Non-mydratic, color image
Field of view	≥ 45 degrees	
Illumination for retina image (Capture)	White LED	Strobe mode with flashing illumination.
Cornea Image (Capture)	White LED	Strobe mode with flashing illumination
Illumination during alignment to patient's retina	NIR LED	Central wavelength in the range of 735-850nm
Focus Diopter adjustment range	-15D to +10 D	Without compensation lens
	-30D to -10D or +5D to +30D	With compensation lens
Minimum pupil size	4 mm	
Focus Adjustment	Auto/ Manual	Split-image technique
Image sensor	CMOS 12 Megapixel	
Z-ranging (Working distance)	2 fiber dots	
Working Distance	25mm from lens to cornea	Accuracy: +/- 0.5mm
Fixation	Internal	10 points

General

Function	Value/Type	Remark
Alignment	Fully automatic 3D tracking	

Alignment Mode	Full Auto / Auto/ Manual	
Chinrest	Motorized	
Interface	USB 2.0 port, Lan, HDMI	
Input/ Output format	Image format: JPEG, PNG , DICOM (optional)	
Display	10.1" LCD monitor, touch panel	
Operation Range	Front / Back: 40mm	
	Left/ Right: 90mm	
	Up/ Down: 30mm	
Chinrest Range	Up/ Down: 70mm	

Note-1:

USB interface is used to connect USB mass storage device.

5.2 Environmental conditions

1) Operating Conditions:

- Temperature: 10°C to 35°C
- Humidity: 30% to 90%RH
- Atmospheric pressure: 800–1060 hPa

2) Storage Conditions:

- Temperature: -10–55°C
- Relative Humidity: 10–95%RH
- Atmospheric pressure: 700–1060 hPa

3) Transport Conditions:

- Temperature: -40–70°C
- Relative Humidity: 10–95%RH
- Vibration, Sinusoidal: 10–500Hz, 0.5g
- Shock: 1/2 Sine Wave, 6 msec,
30G peak (packaged)
- Bump: 1/2 Sine Wave, 6 msec,
10G peak (packaged)

5.3 Electric Rating

- Source voltage :AC100-240V
- Frequency :50-60Hz
- Power input : < 150VA



5.4 Fuse specification:

- Rating: 1.6A/250VAC
- Package: 5 mm x20 mm
- Type: Slow blow

6 Maintenance

6.1 Lens Cleaning

It is recommended to regularly clean the Ocular Lens of the EFC-2600 on weekly basis or when needed.

6.1.1 Material required for Ocular Lens:

- a) Diluted Acetone or Lens cleaning solution
- b) Lens cleaning paper

6.1.2 Method:

Wet the lens paper with cleaning solution and wipe the Ocular Lens with one pass in one direction. Discard the used lens paper. Use a new sheet for each repeat cleaning until the Ocular Lens is clean.

6.2 Chinrest and forehead rest

Soak the cleaning cloth or towel in disinfecting solution or use a wet isopropyl alcohol cleaning paper pad. Wipe the chinrest and forehead rest with the cleaning towels or paper pad before or after use.

If the chinrest paper is used, remove one piece for each patient.

6.3 LCD Monitor/Touch panel

Turn off the power first, and use a soft cleaning cloth to wipe the exterior of the LCD display lightly.

Do not press the LCD monitor using an object with a hard tip. Scratches or failure of the LCD monitor may result.



7 Installation

7.1 Forehead adaptor installation

Step 1. Take the forehead adaptor




Step 2. Put into forehead rest

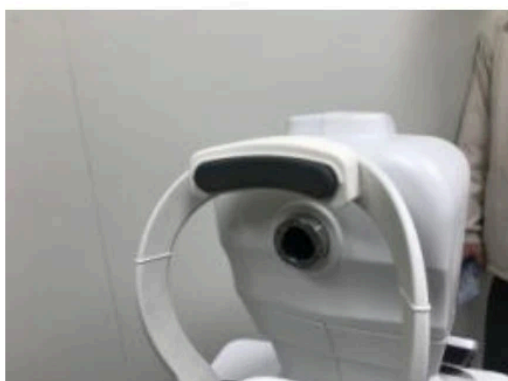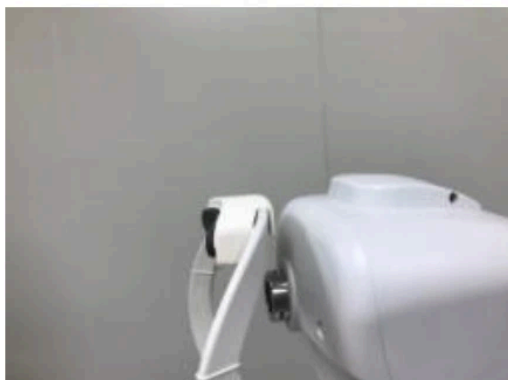


8 Software Update


Step 1. Plug in the USB drive and copy the EFC-2600 installation package file to Windows desktop.

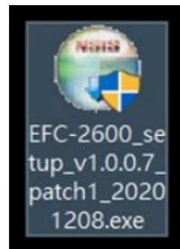


Step 2. Perform the EFC-2600 installation package, below installation dialog will be shown.

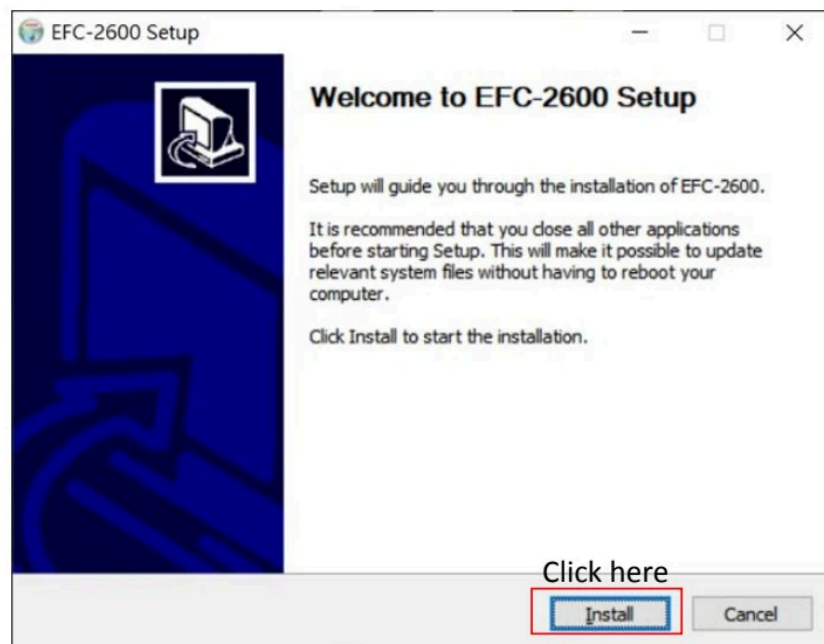

Step 3. Click the “Install” button and wait for installation

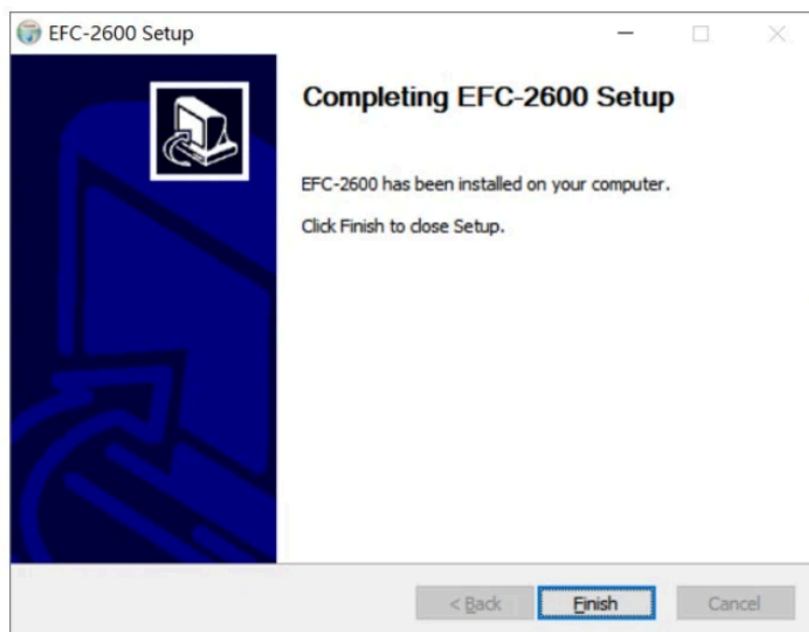

Step 4. Perform the EFC-2600 program.

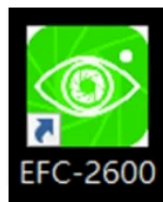



9990 NW 14th Street Unit # 105 Doral, Florida, 33172

United State